


Competitive innovation



Sanat-R 2


Modular revision hip prosthesis

Competitive innovation

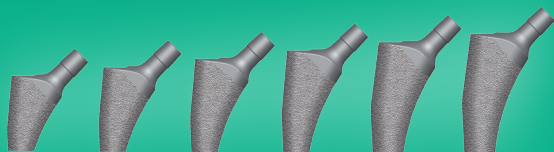

Proximal block

- Wedge shape, 6 different sizes, offset increasing by size
- Titanium raw material, TPS coating for excellent bone ingrowth

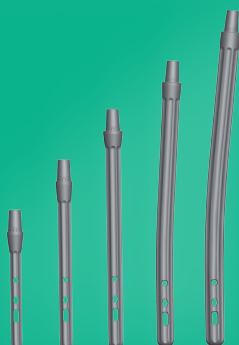

Stem

- Longitudinal ribs for increased stability
- Short-straight and long, anatomically formed stems
- 5 different diameters (12-14-16-18-20 mm)
- Special distal hole for easier distal locking



Unique instruments in endoprothetics

- Fine tuning capability of targeting arm - short stems can be targeted without X-ray
- Hand targeting device for distal targeting of the anatomical stems
- Guide rod assisted rasp introduction



6 different revision flanges for mounted revision shells